

SAM
MEDICAL

SAM ThoraSite®

**POSITION.
ALIGN.
VENT.**

NEW SINGLE SELLING UNIT

- Makes kitting with other products easier
- CE marking and 33 language eIFU provides more locations for use
- Packaged for individual sales




DEVICE TECHNOLOGY

ThoraSite is an anatomical landmark guide, designed to provide a safer approach to the pleural space by pinpointing the procedure site.^{1,2} The device aids in identifying the appropriate intercostal space for lateral thoracostomies and similar procedures. ThoraSite is non-sterile and intended to be single use.


SAFE ZONE WINDOW PRINCIPLE

ThoraSite's patented design references the patient's own anatomical landmarks to position the Safe Zone Window directly over the ideal procedure site. The principal utilizes two axes to create an intersection to safely orient the device. The first axis positions ThoraSite accurately between the anterior axillary and midaxillary line, while the second axis aligns the window over the 3rd, 4th, or 5th intercostal space. This easily guides the user over the appropriate site to vent.

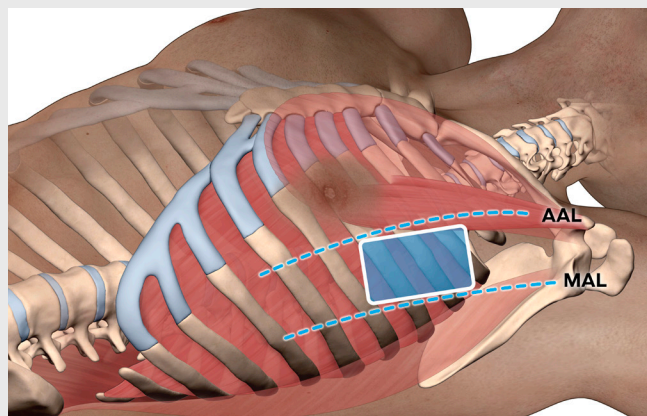

ENGINEERED FOR ACCURACY

- A SAFE ZONE WINDOW**
Highlighted alignment area safely guides user to the appropriate procedure site.
- B AXILLA HOOK**
Developed to be inserted into the patient's axilla for device alignment between the anterior axillary line (AAL) and the midaxillary line (MAL).
- C ALIGNMENT ARROW**
Designed to align the device with the patient's anterior iliac crest/hip.
- D PERFORATED RELEASE**
Created for quick removal when discarding to minimize interference with other devices or procedures.
- E BILATERAL APPLICATION**
Dual sided device accommodates procedures on both patient's left and right side.
- F LOW LIGHT OPTIMIZED**
Printed with high-contrast colors (black/yellow) for better visibility in low-light situations.
- G COMPACT & COMPATIBLE**
Easily fits inside IFAKs and is compatible with all thoracostomy devices.
- H DURABLE FLEX**
Constructed to bend and adapt to patient anatomy while maintaining the Safe Zone Window spatial integrity.

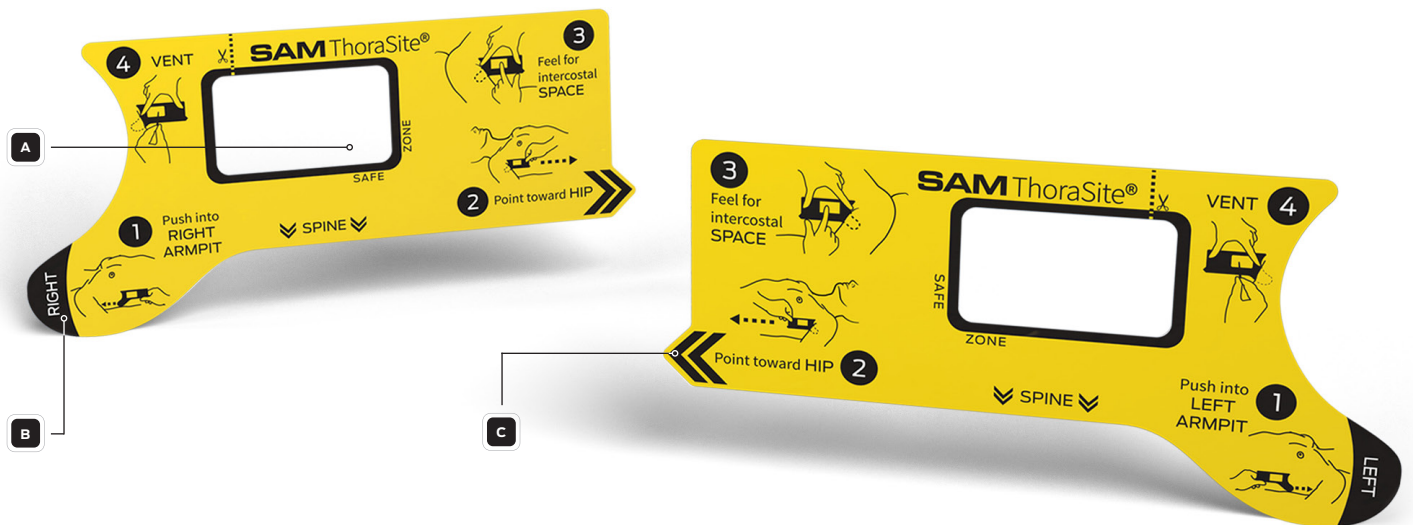

1,2 Shah, Amit Navin MD, MHSc; Kothera, Curt Steven PhD; Dheer, Sachin MD. ThoraSite: A device to improve accuracy of lateral decompression needle and chest tube placement. Journal of Trauma and Acute Care Surgery. 2019;87(1S):S128-S131.

Wernick B, Hon HH, Mubang RN, Cipriano A, Hughes R, Rankin DD, Evans DC, Burfeind WR Jr, Hoey BA, Cipolla J, Galwankar SC, Papadimos TJ, Stawicki SP, Firstenberg MS. Complications of needle thoracostomy: A comprehensive clinical review. Int J Crit Illn Inj Sci. 2015 Jul-Sep;5(3):160-9

PRODUCT SPECIFICATIONS

MODELS AVAILABLE

SAM ThoraSite
TS200-1P
NSN 6515-01-694-8359
While Supplies Last
5-PACK (TS200-5P-EN)
25-PACK (TS200-25P-EN)

PACKAGED SIZE (x1 DEVICE)

Height: 3.0 in (76.20 mm)
Width: 0.04 in (1.02 mm)
Length: 6.125 in (155.60 mm)
Weight: 0.18 oz. (5.1 g)

GLOBAL SALES

sammedical.com
+1 503.639.5474



US Patent N° 10,595,898

SAM® is a registered trademark in the United States and other countries.
©2021 SAM Medical Products. All rights reserved.