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



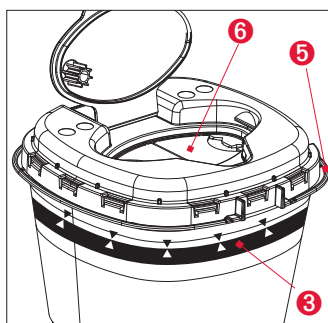
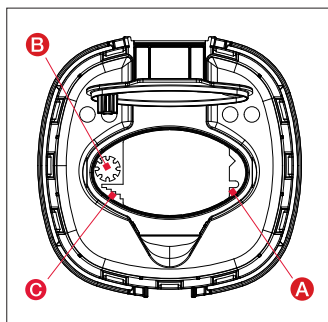
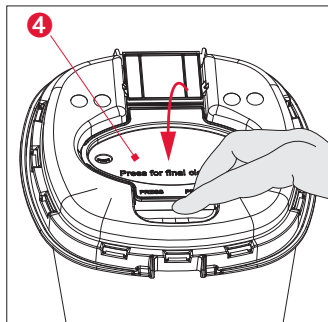
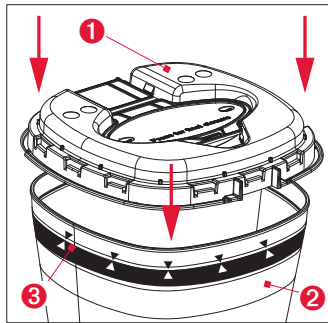
Instructions for use – SARSTEDT Multi-Safe Disposal Boxes Twin Plus

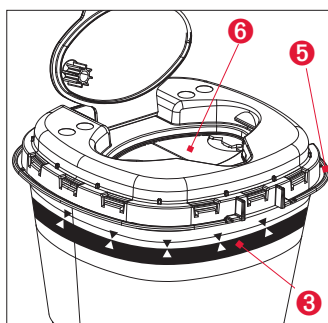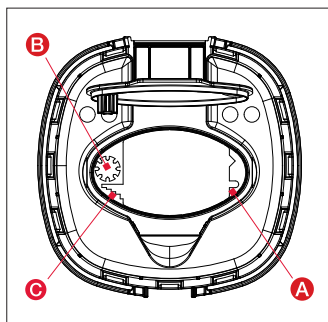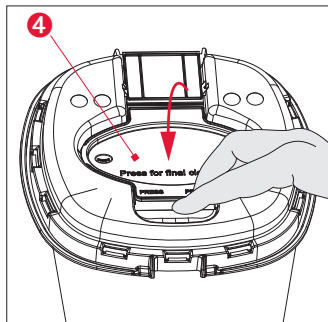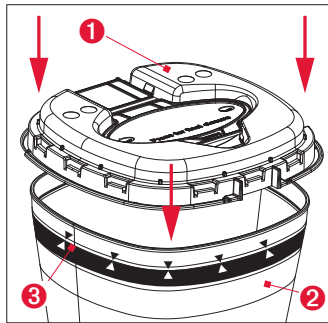1. Intended use

The disposal container is intended exclusively for the disposal of sharp and pointed objects such as cannulas and scalpel blades from the medical sector.

2. Safety instructions:

- Only use boxes of a size suitable for containing the items to be disposed of.
- Keep the box out of the reach of children.
- The lid (1) must be fitted and locked into place before filling begins.
- Attach the box to the recommended adhesive adapter by unscrewing it or secure it by hooking it onto the wall bracket to prevent it from tipping over.
- Only drop items to be disposed of vertically into the box.
- Do not use the daily lid (4) to push the items to be disposed of into the box.
- Scalpels must be disposed of in the box with particular care. Excessive force when throwing items in or refilling with other items poses a risk of the box walls or base becoming warped or damaged.
- Do not force any items into the box.
- Do not pour any liquids into the box.
- Do not fill the box above the fill line (3) (risk of injury!).
- Do not reach into the box with your hand or in any other way (risk of injury!).
- Do not throw the box down.
- Do not shake the box. Do not drop the box.
- Before closing the box, ensure that no objects are protruding from the opening.
- Close the box using only your hands and do not use force.
- Do not press the box against your body (e.g. chest) when closing it.
- Do not squeeze the filled box from the outside.
- Use the carrying handle (5).
- Before disposal, check that the lid (1) and the base (2) are securely closed.

3. Assembly (lid assembly)

Before filling, place the box on a firm surface and press the lid (1) firmly onto the box base (2); the lid (1) will click into place audibly and visibly. The filling opening can be gradually enlarged by breaking out the two segments (6) before fitting the lid.

4. Daily use

For daily use, press the daily lid (4) downwards until it clicks audibly into place.

The lid can be opened by pressing the protruding tab upwards.

5. End cap of the daily lid

After filling the box (Caution: Do not fill above the fill line (3)), proceed as follows: Place the box on a firm surface and press the daily lid (4) firmly into the filling opening until it clicks audibly and visibly all the way round. There is a carrying handle (5) in the upper part of the box, which is bent upwards by pulling gently. The box can now be carried comfortably to the nearest disposal point.

6. Fill line

A fill line (3) is visible on the box label. Do not fill the box beyond this fill line (3).

For safety reasons, the maximum fill level (3) must not be exceeded.

7. Mounting brackets

Adhesive adapters and wall brackets are available to secure the box against impact

8. Twist-off/pull-off mechanism

- 8.1 U-shaped recess (A) for removing Luer systems.
- 8.2 Star-shaped multi-tooth opening (B) as a twist-off mechanism for pens.
- 8.3 Recess with steps (C) for unscrewing cannulas.

9. Disposal and labelling

The disposal boxes are made of environmentally neutral polypropylene. No harmful substances are produced during incineration. Please observe the relevant legal regulations when disposing of the boxes!

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