

EN



TQD

EMS & TACTICAL GEAR



ABOUT THE COMPANY



tqdnipro.com/en

Founded in 2005 as a textile manufacturing company, TQ DNIPRO™ evolved in 2022 into a specialized developer and manufacturer of tactical medical and military equipment in response to the critical need for effective hemorrhage control solutions among military and civilian medics.

Today, TQ DNIPRO™ is a Ukrainian manufacturer of high-quality medical and tactical equipment, combining innovation, reliability and extensive operational experience.

Our products are designed for military personnel, emergency and rescue services, medical professionals and civilian users. Our portfolio includes tourniquets, pouches, first aid kits, stretchers and other specialized equipment designed to meet modern standards and the demands of operational environments.

Our mission

To save lives by developing innovative, combat-proven medical devices and solutions that enable effective casualty care in the most critical conditions.





MANUFACTURED IN UKRAINE

Supporting a Ukrainian manufacturer is a contribution to the country's economy and job creation.



HIGH QUALITY

TQ DNIPRO tourniquets are certified and comply with strict European standards.



INNOVATIVE TECHNOLOGIES

Continuous quality improvement and implementation of innovations to enhance safety levels.



FLEXIBLE COLLABORATION TERMS

We offer special terms, ensuring maximum benefit and convenience for our partners.



QUICK ORDER PROCESSING

We guarantee quick and reliable order processing, with timely dispatch for both domestic and international deliveries.



SUPPORT AT EVERY STAGE

Our clients receive personalized service and continuous support at every stage of cooperation, enabling us to build long-term and mutually beneficial relationships.

TESTING OF TQ DNIPRO TOURNIQUETS HAS SHOWN



ATRAUMATIC

No additional trauma to soft tissues or the neurovascular bundle. Low level of associated pain.



RELIABLE

Maintains the required pressure level and is securely locked against accidental release under load. Operates effectively in extreme climatic conditions.



EFFECTIVE

Ensures complete vascular occlusion and bleeding control. Suitable for different limb sizes, regardless of gender.



FAST

Can be applied to a limb with one hand. Compression is achieved from the first half-turn of the windlass, stopping critical bleeding immediately.

CERTIFICATES AND PATENTS OF TOURNIQUETS TQ DNIPRO GENERATION 2



TQ DNIPRO™ Gen 2 tourniquets are protected by patents:

Ukraine: No. 151561, No. 152755

European Union: No. EP4514238B1, No. EP4456803B1


United States: No. US12349921B2



The manufacturer holds international recognition and certification:

European Union: Declaration of Conformity and documentation in accordance with Regulation (EU) 2017/745. Both the manufacturer and the product are registered in EUDAMED (SRN: UA-MF-000039898).

Australia: Included in the Australian Register of Therapeutic Goods (ARTG) in accordance with the Therapeutic Goods Act 1989 and Therapeutic Goods (Medical Devices) Regulations 2002. United Kingdom: Registered with the MHRA in compliance with the UK Medical Devices Regulations 2002, based on EU technical documentation and the EU Declaration of Conformity.

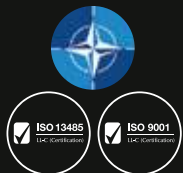

Additional compliance highlights:

Assigned NATO Stock Number (NSN) 6515-61-017-6778

Manufacturer holds an NCAGE code within the NATO codification system

Certified quality management systems: ISO 13485:2016 (Medical Devices) and ISO 9001:2015 (Quality Management System)



**TOURNIQUET TQ DNIPRO GEN 2, BLACK
TACTICAL**

CODE: TQD01

Intended purpose - providing first aid to stop critical
bleeding from the upper and lower limbs:

- during military operations
- in extreme and emergency situations
- in camping conditions
- in ambulances
- in medical institutions
- in everyday life



video review



TOURNIQUET TQ DNIPRO GEN 2, BLUE TRAINING

CODE: TQD02

Intended purpose -the training tourniquet is intended to practice vital skills for stopping critical bleeding from the upper and lower limbs:



- to save your own life and the life of another person
- in emergency situations
- In ambulances, in medical facilities
- In everyday life



TOURNIQUET TQ DNIPRO GEN 2, ORANGE
EMERGENCY MEDICAL SERVICES

CODE: TQD03

Intended use – providing first aid to stop critical bleeding of the upper and lower extremities in emergency situations:

- to save your own life and the life of another person
- in emergency situations
- In ambulances, in medical facilities
- In everyday life



The image displays three tourniquet pouches against a dark background. The top pouch features a camouflage pattern and has a small olive green label with 'TG ONPRD' printed on it. The middle pouch is olive green and also has a small olive green label with 'TG ONPRD'. The bottom pouch is bright red and has a small black label with 'TG ONPRD' in white text. A central text box with a gold border contains the text 'TOURNIQUET POUCH MODEL1 GEN 2'.

TOURNIQUET POUCH
MODEL1 GEN 2

TG ONPRD

VENDOR CODE:
PM019 GEN2



COLOR: COYOTE

VENDOR CODE:
PM011 GEN2



COLOR: OLIVE

VENDOR CODE:
PM012 GEN2



COLOR: BLACK

VENDOR CODE:
PM0112-GEN2



COLOR: NAVY

VENDOR CODE:
PM018 GEN2



COLOR: MULTICAM

VENDOR CODE:
PM0110 GEN2



COLOR: RED

DESIGN FEATURES

DURABILITY

Made of lightweight and durable materials that protect the tourniquet from mechanical damage.

FUNCTIONALITY

Versatility. Side walls made of elastic straps allow storing not only tourniquets but also other similarly sized items, making the holder a universal accessory. Attaches to equipment and provides additional accessory fixation with elastic cords on the sides.

ATTACHMENT METHOD:

The holder is attached to the gear using the MOLLE fastener

TECHNICAL SPECIFICATIONS

MATERIAL:

Elastic strap – Italy
Polyamide sling – Ukraine
(USA for Multicam color)

Weight: 30 g

Size: 175 x 60 x 25 mm

TOURNIQUET POUCH
MODEL 2

VENDOR CODE:
PM025



COLOR: COYOTE

VENDOR CODE:
PM021



COLOR: OLIVE

VENDOR CODE:
PM022



COLOR: BLACK

VENDOR CODE:
PM028



COLOR: MULTICAM



TECHNICAL SPECIFICATIONS

MATERIAL:
CORDURA® Squadron® laminate
Germany
(USA for Multicam color)
Elastic strap – Italy

Weight: 80 g
Size: 225 x 50 x 40 mm

DESIGN FEATURES

PROTECTION

Fully enclosed flap-type design provides reliable protection of the tourniquet from environmental exposure, ultraviolet radiation, moisture, dirt, small debris, and mechanical damage during transport and field operations.

FUNCTIONALITY

Designed for storage and transportation of windlass-type tourniquets with C-shaped retainers. Equipped with dedicated attachment points for tactical shears, a permanent marker for recording tourniquet application time, and examination gloves.

ATTACHMENT METHOD:

Features a universal MOLLE mounting system compatible with military, law enforcement, and civilian equipment. An additional fixation system allows attachment to belts and various load-bearing platforms.

TOURNIQUET POUCH
PM03

VENDOR CODE:
PM035



COLOR: COYOTE

VENDOR CODE:
PM031



COLOR: OLIVE

VENDOR CODE:
PM032



COLOR: BLACK

VENDOR CODE:
PM0312



COLOR: NAVY

VENDOR CODE:
PM038



COLOR: MULTICAM

VENDOR CODE:
PM0310



COLOR: RED

DESIGN FEATURES

DURABILITY

Made of lightweight and durable materials that protect the tourniquet from mechanical damage.

FUNCTIONALITY

Versatility. Side walls made of elastic straps allow storing not only tourniquets but also other similarly sized items, making the holder a universal accessory. Attaches to equipment and provides additional accessory fixation with elastic cords on the sides.

ATTACHMENT METHOD:

The holder is attached to the gear using the MOLLE fastener

TECHNICAL SPECIFICATIONS

MATERIAL:

Cordura fabric - Germany
(USA for Multicam color)
Elastic strap - Italy

Weight: 30 g

Size: 175 x 60 x 25 mm

TOURNIQUET POUCH
PM08 GEN2



VENDOR CODE:
PM085 GEN2



COLOR: COYOTE

VENDOR CODE:
PM081 GEN2



COLOR: OLIVE

VENDOR CODE:
PM082 GEN2



COLOR: BLACK

VENDOR CODE:
PM0812-GEN2



COLOR: NAVY

VENDOR CODE:
PM088 GEN2



COLOR: MULTICAM

VENDOR CODE:
PM0810 GEN2



COLOR: RED

DESIGN FEATURES

DURABILITY

Made from lightweight and durable materials that effectively protect the tourniquet from UV rays, moisture, and mechanical damage.

FUNCTIONALITY

The internal intermediate latch securely fixes the tourniquet from falling out during opening. It is easy to open if necessary to remove the tourniquet, which reduces the time of applying the tourniquet separating from the "body" of the tourniquet the part with which the hand comes into contact.

The quick opening system of the side flaps together with the valve, which opens completely, provides instant access to the turnstile.

Free positioning of the tourniquet with a C-shaped latch.

ATTACHMENT METHOD:

Equipped with a MOLLE mounting system for attaching tourniquets to military, law enforcement, or civilian gear. An additional fixation system allows attachment to equipment without MOLLE compatibility.

TECHNICAL SPECIFICATIONS

MATERIAL:

Cordura fabric – Germany
(USA for Multicam color)
Elastic strap – Italy

Weight: 90 g
Size: 230x50x40 mm



A close-up photograph of a soldier's leg in camouflage pants. A tan-colored IFAK (Individual First Aid Kit) bag is attached to the leg using a brown webbing strap with a metal buckle. The bag has a smaller front pocket and a Velcro patch. Above the bag, another strap with three metal rings is visible. A semi-transparent dark green box with white text is overlaid on the left side of the image.

IFAK BAG
PM 21 GEN2



VENDOR CODE:
PM215 GEN2

COLOR: COYOTE



VENDOR CODE:
PM211 GEN2

COLOR: OLIVE



VENDOR CODE:
PM212 GEN2

COLOR: BLACK



VENDOR CODE:
PM2112-GEN2

COLOR: NAVY



VENDOR CODE:
PM218 GEN2

COLOR: MULTICAM



VENDOR CODE:
PM2110 GEN2

COLOR: RED

DESIGN FEATURES

DURABILITY

Lightweight and durable materials effectively protect the contents from moisture, UV rays, and mechanical damage.

ATTACHMENT METHOD

The pouch attaches via the MOLLE system to the platform and is secured to the thigh with elastic straps and a FAST buckle.

FUNCTIONALITY

Includes an organizer, a drainage hole for excess moisture, and an external pocket. Ensures a quick and easy opening with access to all medical supplies without detaching from the platform.

TECHNICAL SPECIFICATIONS

MATERIAL:

Cordura fabric – Germany

(USA for Multicam color)

Elastic strap – Italy.

FAST polyamide buckle – Italy

Weight: 510 g

Size: 145 × 180 × 80 mm



IFAK BAG
PM 22 GEN2

VENDOR CODE:
PM225 GEN2



COLOR: COYOTE

VENDOR CODE:
PM221 GEN2



COLOR: OLIVE

VENDOR CODE:
PM222 GEN2



COLOR: BLACK

VENDOR CODE:
PPM2212-GEN 2



COLOR: NAVY

VENDOR CODE:
PM228 GEN2



COLOR: MULTICAM

VENDOR CODE:
PM2210 GEN2



COLOR: RED

DESIGN FEATURES

DURABILITY

Lightweight and durable materials effectively protect the contents from moisture, UV rays, and mechanical damage.

ATTACHMENT METHOD

The pouch attaches to the platform via Velcro.
The platform with the pouch is secured to the gear using the MOLLE system.

FUNCTIONALITY

Features a convenient internal organizer and a drainage hole at the bottom for excess liquid and moisture. The design ensures quick opening and easy access to all medical supplies without detaching the pouch from the platform. Can be easily detached from the base for use or relocation.

TECHNICAL SPECIFICATIONS

MATERIALS:

Cordura fabric – Germany

(USA for Multicam color)

Elastic strap – Italy

FAST polyamide buckle – Italy

Weight: 310 g

Size: 145 × 180 × 80 mm



IFAK BAG
PM 23

VENDOR CODE:
PM235



COLOR: COYOTE

VENDOR CODE:
PM231



COLOR: OLIVE

VENDOR CODE:
PM232



COLOR: BLACK

VENDOR CODE:
PM2312



COLOR: NAVY

VENDOR CODE:
PM238



COLOR: MULTICAM

VENDOR CODE:
PM2310



COLOR: RED

DESIGN FEATURES

DURABILITY

Lightweight and durable materials protect the contents from UV rays and mechanical damage.

ATTACHMENT METHOD

Designed for secure retention on the lower leg (shin).
The pouch is attached to the shin using Velcro and secured with a FAST buckle.

FUNCTIONALITY

Provides quick and easy access to all medical supplies without detaching the pouch from the shin.

TECHNICAL SPECIFICATIONS

MATERIAL:

Cordura fabric – Germany
(USA for Multicam color)
Elastic strap – Italy
FAST polyamide buckle – Italy

Weight: 135 g

Size: 510 × 120 × 40 mm



CLASSIFICATION:

Full IFAK (Individual First Aid Kit)

Designed for shin-mounted carry, allowing the user to keep essential medical supplies accessible without overloading the upper-body equipment. The pouch provides quick access to the contents while maintaining freedom of movement and can be rapidly removed when necessary.



IFAK BAG
PM 24 GEN2

VENDOR CODE:
PM245 GEN2



COLOR: COYOTE

VENDOR CODE:
PM241 GEN2



COLOR: OLIVE

VENDOR CODE:
PM242 GEN2



COLOR: BLACK

VENDOR CODE:
PM2412-GEN 2



COLOR: NAVY

VENDOR CODE:
PM248 GEN2



COLOR: MULTICAM

VENDOR CODE:
PM2410 GEN2



COLOR: RED

DESIGN FEATURES

DURABILITY

Lightweight and durable materials effectively protect the contents from moisture, UV rays, and mechanical damage.

ATTACHMENT METHOD

The pouch attaches to the gear using the MOLLE system, which allows for quick removal.

FUNCTIONALITY

Features a convenient internal organizer for medical supplies and a drainage hole at the bottom for excess liquid and moisture.
The design ensures quick opening and easy access to all medical supplies without detaching the pouch from the platform.
Can be quickly and easily removed from the gear for use or relocation.

TECHNICAL SPECIFICATIONS

MATERIALS:

Cordura fabric – Germany
(USA for Multicam color)
Elastic strap – Italy
FAST polyamide buckle – Italy

Weight: 310 g
Size: 150 × 200 × 90 mm

IFAK BAG
PM 25 GEN2





VENDOR CODE:
PM255 GEN2

COLOR: COYOTE



VENDOR CODE:
PM251 GEN2

COLOR: OLIVE



VENDOR CODE:
PM252 GEN2

COLOR: BLACK



VENDOR CODE:
PM2512-GEN2

COLOR: NAVY



VENDOR CODE:
PM258 GEN2

COLOR: MULTICAM



VENDOR CODE:
PM2510 GEN2

COLOR: RED

DESIGN FEATURES

DURABILITY

Lightweight and durable materials effectively protect the contents from moisture, UV rays, and mechanical damage.

ATTACHMENT METHOD

The pouch can be secured to the thigh with adjustable height settings. The attachment system includes a vertical FAST buckle for securing to a belt and a horizontal elastic strap with a FAST buckle to prevent lateral movement.

FUNCTIONALITY

The optimized internal layout allows for organizing medical supplies according to echelon principles, ensuring quick access in critical situations.

TECHNICAL SPECIFICATIONS

MATERIAL:

Cordura fabric – Germany
(USA for Multicam color)
Elastic strap – Italy
FAST polyamide buckle – Italy

Weight: 410 g

Size: 410 × 120 × 40 mm





IFAK BAG
PM 26

VENDOR CODE:
PM265



COLOR: COYOTE

VENDOR CODE:
PM261



COLOR: OLIVE

VENDOR CODE:
PM262



COLOR: BLACK

VENDOR CODE:
PPM2612



COLOR: NAVY

VENDOR CODE:
PM268



COLOR: MULTICAM

VENDOR CODE:
PM2610



COLOR: RED

DESIGN FEATURES

RAPID ACCESS

The tear-away platform enables immediate deployment of the IFAK when every second matters. The pouch can be detached with a single motion for rapid self-aid or treatment of others.

FUNCTIONALITY

The dual-zip opening provides wide access to medical supplies, while the organized interior layout keeps essential equipment secure, visible, and quickly accessible under pressure.

ATTACHMENT METHOD

The pouch is mounted to gear using a MOLLE-compatible platform. The hook and loop fixation ensures reliable retention during movement while allowing fast removal when needed.

TECHNICAL SPECIFICATIONS

MATERIALS:

Cordura fabric – Germany
(USA for Multicam color)

Elastic strap – Italy

FAST polyamide buckle – Italy

Weight: 310 g

Size: 145 × 180 × 80 mm

IFAK POUCH
PM 27 GEN 2



VENDOR CODE:
PM275 GEN2



COLOR: COYOTE

VENDOR CODE:
PM271 GEN2



COLOR: OLIVE

VENDOR CODE:
PM272 GEN2



COLOR: BLACK

VENDOR CODE:
PM2712-GEN2



COLOR: NAVY

VENDOR CODE:
PM278 GEN2



COLOR: MULTICAM

VENDOR CODE:
PM2710 GEN2



COLOR: RED

DESIGN FEATURES

DURABILITY

Lightweight and durable materials protect the contents from moisture, UV rays, and mechanical damage.

ATTACHMENT METHOD

Can also be attached to gear using the MOLLE system.

FUNCTIONALITY

Can be quickly extracted from the case from either side.
Equipped with two external mounts for tourniquets.

TECHNICAL SPECIFICATIONS

MATERIALS:

Cordura fabric – Germany

(USA for Multicam color)

Elastic strap – Italy

FAST polyamide buckle – Italy

Weight: 240 g

Size: 220 × 190 × 60mm



CLASSIFICATION:

Full IFAK (Individual First Aid Kit)

The outer pouch allows the placement of two tourniquets (top and bottom of the pouch) for quick and easy access.

The direction of the tourniquets can be each for a different hand.

IFAK POUCH
PM 28

The image shows a person wearing a camouflage tactical vest. A large, rectangular IFAK pouch is attached to the side of the vest. The pouch is made of camouflage fabric and features two large, silver-colored buckles. A black strap is visible at the top of the pouch. The vest itself has multiple pockets and straps, some of which are labeled "MUTICAM". The person is wearing a tan-colored shirt and camouflage pants.

VENDOR CODE:
PM285



COLOR: COYOTE

VENDOR CODE:
PM281



COLOR: OLIVE

VENDOR CODE:
PM282



COLOR: BLACK

VENDOR CODE:
PM2812



COLOR: NAVY

VENDOR CODE:
PM288



COLOR: MULTICAM

VENDOR CODE:
PM2810



COLOR: RED

DESIGN FEATURES

RAPID DEPLOYMENT

Tunnel-type construction with a removable internal organizer. The organizer can be extracted from either side using the left or right hand, ensuring fast access to medical supplies in emergency situations.

FUNCTIONALITY

The outer sleeve includes holders for two windlass tourniquets, a dedicated trauma shears pocket, retractor attachment loops, and a hook and loop panel for identification patches. The internal organizer allows structured placement of medical supplies according to the MARCH protocol.

ATTACHMENT METHOD

Compatible with MOLLE equipment. The removable platform enables attachment both to tactical gear and to a standard belt for versatile carry options.

TECHNICAL SPECIFICATIONS

MATERIALS:

Cordura fabric – Germany

(USA for Multicam color)

Elastic strap – Italy

FAST polyamide buckle – Italy

Weight: 390 g

Size: 190 × 180 × 80mm

CLASSIFICATION:

Full IFAK (Individual First Aid Kit)


The tunnel-style pouch features a removable medical insert that can be deployed from either side, allowing rapid access with either hand.

The outer pouch includes dedicated holders for two windlass tourniquets and a trauma shears pocket, ensuring immediate access to critical life-saving equipment.

The internal organizer is designed for structured placement of medical supplies according to the MARCH protocol, improving efficiency during emergency response.



IFAK BAG
PM 29

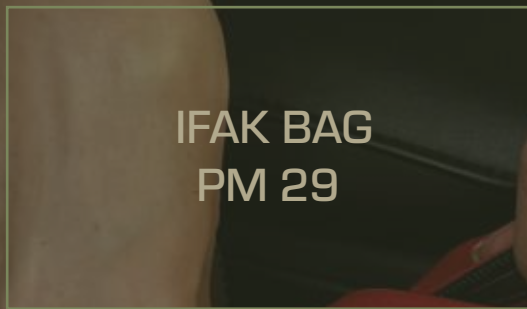

A pouch designed for home, office, or vehicle use. Intended for transporting and storing first aid supplies.

VENDOR CODE:
PM 2910

COLOR: RED

DESIGN FEATURES

DURABILITY

Lightweight and durable materials protect the contents from moisture, UV rays, and mechanical damage.

ATTACHMENT METHOD

Equipped with two handles for easy transportation and opening. Features an external pocket for gloves and outer compartments for scissors, a marker, and a tourniquet.

FUNCTIONALITY

Easy to open.
Provides quick access to medical supplies in critical situations.

TECHNICAL SPECIFICATIONS

MATERIAL:

Cordura fabric – Germany
Elastic strap – Italy
FAST polyamide buckle – Italy

Weight: 385 g
Size: 190 × 145 × 95 mm

A close-up photograph of a soldier wearing a camouflage uniform. The soldier is carrying a large, modular tactical bag, identified as a TQD BM30 medic bag, which is also in camouflage. The bag features multiple straps, buckles, and a mesh pocket. A black bag is attached to the bottom of the main bag. The background is a blurred, dry, brownish landscape.

TQD BM30
MEDIC BAG

Medic Bag TQD30 is a compact, versatile medical bag for casualty care according to the MARCH algorithm.

- Rapid access to critical medical equipment
- Organized central storage compartment
- Customizable configuration
- Soft stretcher attachment and quick deployment system
- Designed for military, civilian, and rescue medics

VENDOR CODE:
BM305



COLOR: COYOTE

VENDOR CODE:
BM301



COLOR: OLIVE



VENDOR CODE:
BM302



COLOR: BLACK

VENDOR CODE:
BM308



COLOR: MULTICAM

TECHNICAL SPECIFICATIONS

Fabric: Cordura® 1000D

(Germany for Olive, Black colors,
USA for Multicam®, Coyote colors)

Elastic bands and mesh – Italy

Polyamide webbing – Ukraine

Hardware – Italy, 2M

Zippers – YKK, Japan Original

Weight: 1,410 g ± 0.02 kg

Volume: 10 L

Dimensions: 400x250x150 mm

TQD BM33
MEDIC BAG
HYBRID

Compact hybrid bag for MARCH-based pre-hospital care, supporting multiple casualties and infusion therapy.

Features:

- Rapid access to critical medical equipment
- Organized storage with external pockets and internal organizers
- Suitable for military and civilian medics, and rescue personnel
- Functions as both a tactical medic bag and an EDC bag.

VENDOR CODE:
BM335



COLOR: COYOTE

VENDOR CODE:
BM331



COLOR: OLIVE



VENDOR CODE:
BM332



COLOR: BLACK

VENDOR CODE:
BM338



COLOR: MULTICAM

TECHNICAL SPECIFICATIONS

Fabric: Cordura® 1000D
(Germany for Olive, Black colors,
USA for Multicam®, Coyote colors)
Elastic bands and mesh – Italy
Polyamide webbing – Ukraine
Hardware – Italy, 2M
Zippers – YKK, Japan Original
Weight: 1,150 g ± 0.02 kg
Volume: 6 L
Dimensions: 330 × 170 × 110 mm



UV
protection



Flame
resistance



Water-
repellent



Durability



Inside
organizer



Ampoule
organizer

A person wearing a black tactical vest and a large black medic bag around their waist. The person is also wearing black gloves and a watch. The background is a dark, industrial setting with concrete walls and metal structures.

TQD BM34
MEDIC BAG

The TQD34 Medic Bag is a compact solution for carrying medical supplies and providing first aid to multiple casualties in accordance with the MARCH algorithm. Designed for military and civilian medics, as well as rescue personnel.

External pockets and internal organizers provide rapid access to hemorrhage control equipment, lighting tools, and casualty assessment devices. The spacious main compartment ensures efficient organization of additional medical supplies.

Its modular design adapts to mission requirements, while the integrated bottom carrying system allows rapid deployment of soft stretchers.

VENDOR CODE:
BM345



COLOR: COYOTE

VENDOR CODE:
BM341



COLOR: OLIVE

VENDOR CODE:
BM342



COLOR: BLACK

VENDOR CODE:
BM348



COLOR: MULTICAM



TECHNICAL SPECIFICATIONS

Fabric: Cordura® 1000D
(Germany for Olive, Black colors,
USA for Multicam®, Coyote colors)
Elastic bands and mesh – Italy
Polyamide webbing – Ukraine
Hardware – Italy, 2M
Zippers – YKK, Japan Original
Weight: 1,410 g ± 0.02 kg
Volume: 10 L
Dimensions: 400x250x150 mm

TQD BM35
MEDIC BAG
HYBRID

Compact medical bag for MARCH-based pre-hospital care.

- Multiple casualty support
- Infusion therapy capable
- Quick-access pocket system
- Organized internal layout
- Tactical and everyday-use design

For military and civilian medics, as well as rescue personnel.



VENDOR CODE:
BM355



COLOR: COYOTE

VENDOR CODE:
BM351



COLOR: OLIVE

VENDOR CODE:
BM352



COLOR: BLACK

VENDOR CODE:
BM358



COLOR: MULTICAM

TECHNICAL SPECIFICATIONS

Fabric: Cordura® 1000D

(Germany for Olive, Black colors,
USA for Multicam®, Coyote colors)

Elastic bands and mesh – Italy

Polyamide webbing – Ukraine

Hardware – Italy, 2M

Zippers – YKK, Japan Original

Weight: 1,150 g ± 0.02 kg

Volume: 6 L

Dimensions: 330 × 170 × 110 mm



UV
protection



Flame
resistance



Water-
repellent



Durability



Inside
organizer



Ampoule
organizer

WATCH PROTECTOR
WP 01

VENDOR CODE:
WP015

COLOR:
COYOTE



VENDOR CODE:
WP012

COLOR:
BLACK



VENDOR CODE:
WP011

COLOR:
OLIVE



VENDOR CODE:
WP018

COLOR:
MULTICAM



VENDOR CODE:
WP0110

COLOR:
RED




PURPOSE



SCRATCHES

Watch Protector is intended to prevent scratches on the dial, however dirt and debris can be accumulated in the material and scratch the surface of the watch and dial.

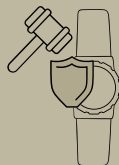

HITS

Even in peaceful life, you may not notice and accidentally break your watch, what to speak of the situations at the front. To prevent such troubles, we recommend using our Tactical Watch Protector.



LIGHT MASKING AND GLARE

If you use a smart watch, then you just need light masking, both for the watch screen and illumination of the sensors. Anti-glare protection is necessary for mechanics and quartz. Our Watch Protector solves both problems.

DESIGN FEATURES

DURABILITY

Made from lightweight and durable materials, 100% polyamide.

ATTACHMENT METHOD

Worn over the watch and secured with a Velcro fastener.

FUNCTIONALITY

Conceals the watch`s screen illumination and sensor lights.
Protects the watch face from glare and potential damage.


TECHNICAL SPECIFICATIONS

MATERIALS:

Cordura fabric 500D – Germany
(USA for Multicam color)
eco-leather

DIMENSIONS:

Strap dimensions:
Length - 275 mm
Width (dial) - 50 mm
Width (strap) - 40 mm
Recommended size
for watch
body - dia 30 - 55 mm

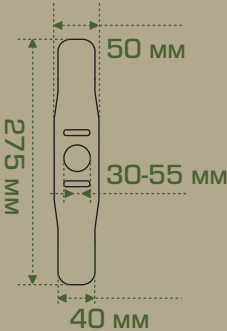



EVACUATION STRETCHERS

VENDOR CODE: STE027



COLOR: SAND

SOFT STRETCHER E

DESIGN FEATURES

DURABILITY

Made of water-repellent, durable fabric. The continuous strap running under the stretchers ensures structural strength. All fabric and strap joints are reinforced with strong stitching for additional durability.

ATTACHMENT

Separate leg and chest fastening ensures body stability during transport, minimizing additional injuries from movement and helps evenly distribute the load to prevent shifting during carrying.

FUNCTIONALITY

When unfolded, the stretchers create a flat horizontal surface for convenient transport. Integrated handles maintain the center of gravity optimally, minimizing the risk of tipping and improving stability during movement.

Handles allow two or more rescuers to coordinate efforts, maintaining balance and control over the injured person.

Handles near the center of gravity reduce physical strain on rescuers' backs and arms.

Proper center of gravity placement makes these stretchers easier to maneuver in confined spaces, enhancing evacuation efficiency.

TECHNICAL SPECIFICATIONS

MATERIALS:

Cordura 1000D fabric
Polyamide strap - Ukraine
Fittings - Italy

DIMENSIONS:

Weight: Up to 1.4 kg
Size in Transport Position:
250 x 220 mm
Size: 2500 x 900 mm

Maximum Load: 140 kg
Number of Handles: 9 handles:
4 on each side along the length,
1 handle at the head end,
4 wrist straps.

SOFT STRETCHER BAG "DNIPRO" MODEL NO. 5

VENDOR CODE: PM057
COLOR: SAND



VENDOR CODE: PM052
COLOR: BLACK



MATERIALS:

Polyamide 100% (water-repellent)
Velcro, Type of fastening: MOLLE

DIMENSIONS:

280 x 210 x 100 mm

SOFT STRETCHER S

DESIGN FEATURES

DURABILITY

Made of water-repellent, durable fabric. The continuous strap running under the stretchers ensures structural strength. All fabric and strap joints are reinforced with strong stitching for additional durability.

ATTACHMENT

Separate leg and chest fastening ensures body stability during transport, minimizing additional injuries from movement and helps evenly distribute the load to prevent shifting during carrying.

FUNCTIONALITY

When unfolded, the stretchers create a flat horizontal surface for convenient transport.

Handles positioned away from the structure provide better maneuverability, essential in confined spaces or when frequent hand changes are needed.

These stretchers are suitable for short-distance transport, particularly in challenging terrains.

TECHNICAL SPECIFICATIONS

MATERIALS:

Oxford 600 fabric

Polyamide strap – Ukraine

Fittings – Italy

DIMENSIONS:

Weight: Up to 1.1 kg

Size in Transport Position: 250 x 220 mm

Size: 2500 x 900 mm

Maximum Load: 140 kg

Number of Handles: 9 handles:

4 on each side along the length,

1 handle at the head end,

4 wrist straps.



VENDOR CODE:
STS022

COLOR: BLACK



VENDOR CODE:
STS027

COLOR: SAND

SOFT STRETCHER L

DESIGN FEATURES

The DNIPRO-L evacuation stretcher is designed for rapid casualty extraction and transportation in tactical, emergency, and disaster-response environments.

Manufactured from ultra-lightweight PTFE-coated fabric, the stretcher combines exceptional strength with minimal weight, allowing rescuers to reduce carried load while maintaining reliable casualty handling capabilities.

Its compact, foldable design enables convenient storage in medical kits, vehicles, and emergency equipment packs. Multiple reinforced carry handles provide efficient casualty movement by one or more rescuers, while integrated retention straps help secure the casualty during evacuation.

The DNIPRO-L evacuation stretcher is suitable for military, law enforcement, rescue, humanitarian, and civilian emergency applications where lightweight, portable evacuation equipment is essential.

VENDOR CODE: STLO25
COLOR: COYOTE



TECHNICAL SPECIFICATIONS

MATERIALS:

Fabric: Teflon
Polyamide strap - Ukraine
Fittings - Italy

DIMENSIONS:

Weight: Up to 0,75 kg
Size: 1900 x 720 mm

Maximum Load: 140 kg
Number of Handles: 9 handles:
4 on each side along the length,
1 handle at the head end,
4 wrist straps.





 LinkedIn: tqdnipro
 +380 67 115 7047
 globalsales@dnipro-tmt.com
 Dnipro, Ukraine



tqdnipro.com/en

